

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Pentoxyverine citrate (Antitussive)

ペントキシベリン クエン酸塩 (鎮咳薬)

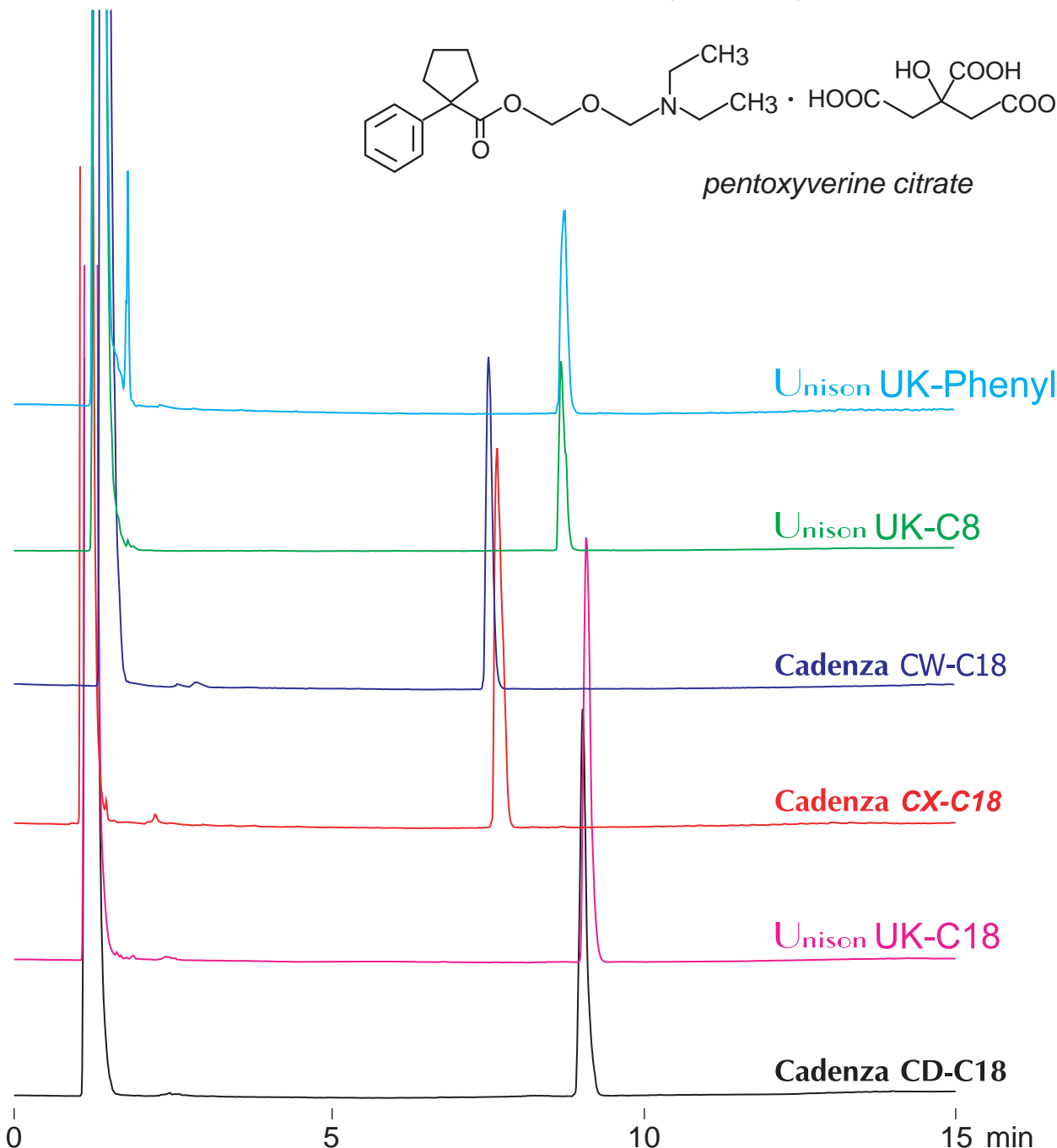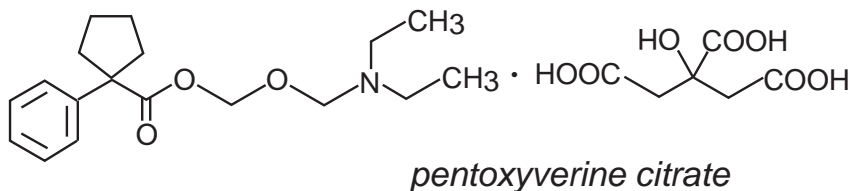

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate= 90 / 10

25-80%B (0-15min)

1.0mL/min (8MPa), 37°C

ELSD (spray chamber 30°C, drift chamber 52°C)

10µL (1mg/mL, methanol)