


Intrada Organic Acid

150 x 2 mm


HPLC COLUMN APPLICATION

LC-MS: Bialaphos, Glyphosate


ビアラホス, グリホサートの LC-MS 分析

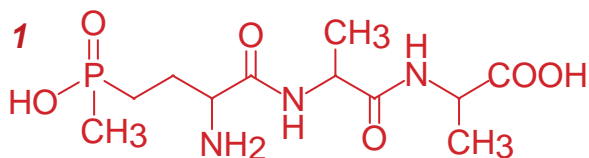

bialaphos

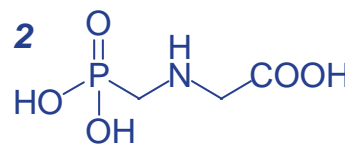

glyphosate

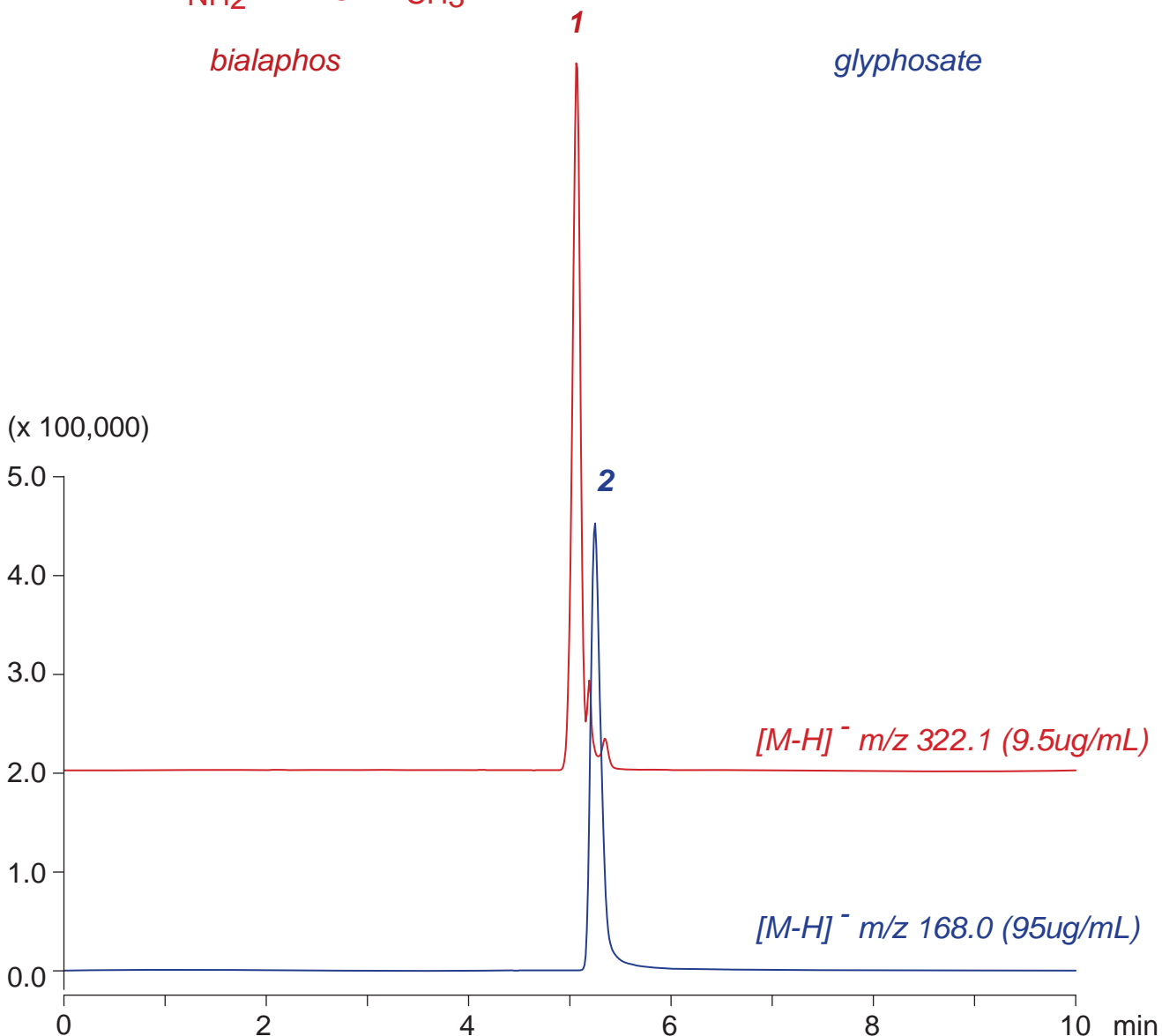

Intrada Organic Acid, 150 x 2 mm

A: acetonitrile / 100mM ammonium formate = 95 / 5

B: acetonitrile / 100mM ammonium formate = 10 / 90

10-100%B (0-7min), 100%B (7-10min)

0.2mL/min (6MPa), 37°C, MS (ESI, SIM, negative)

5µL (9.5-95µg/mL, water / HCOOH = 100 / 2)