


Cadenza **CX-C18**  
 Cadenza **CL-C18** Cadenza **CW-C18**  
 Cadenza **CD-C18** Unison **UK-C18**

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Benzpyrene isomers

ベンゾピレン異性体

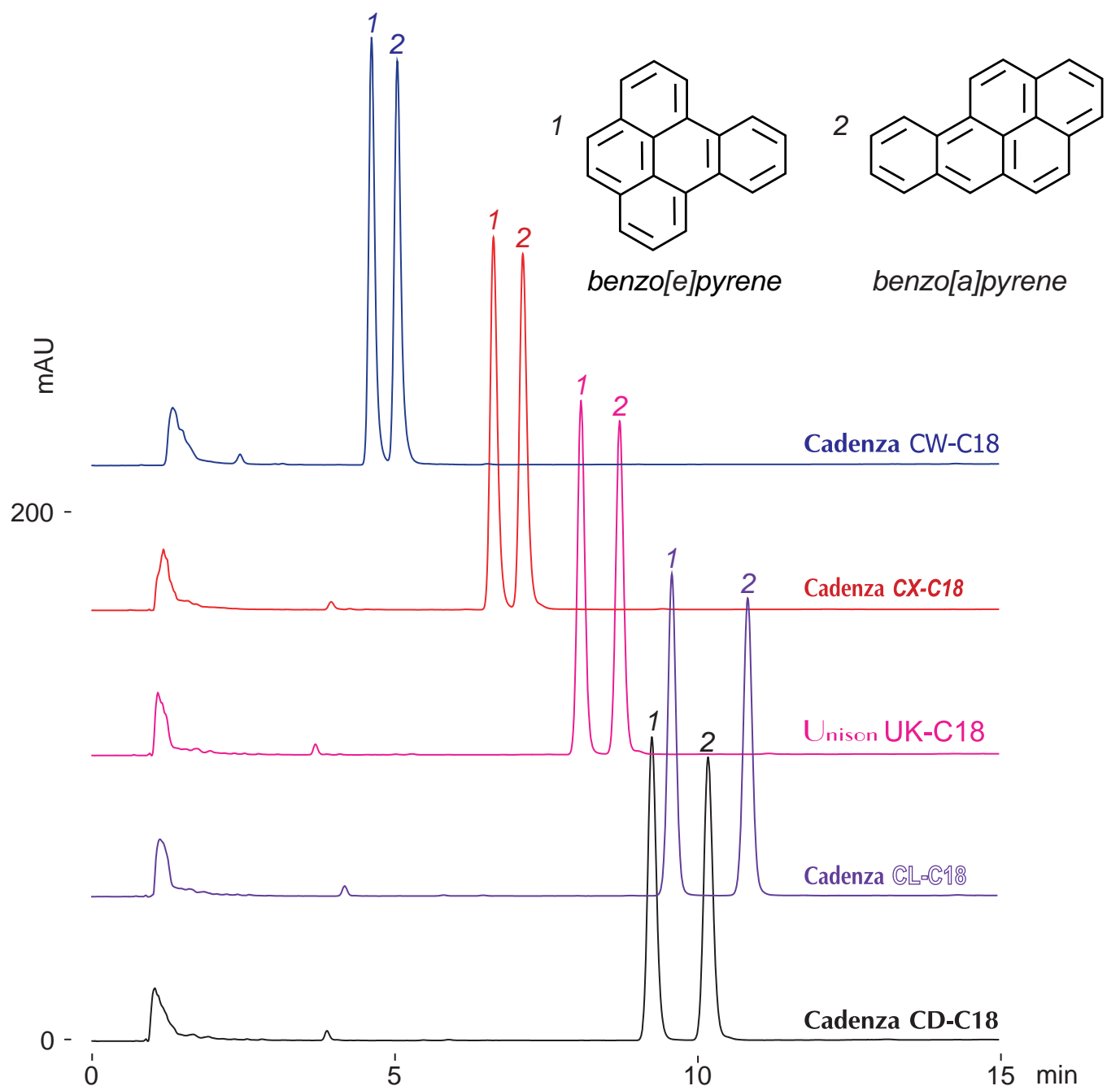

100 x 4.6 mm

A: water / acetonitrile / acetic acid = 50 / 50 / 0.1

B: water / acetonitrile / acetic acid = 10 / 90 / 0.1

65-100%B (0-15min)

1.0mL/min (6MPa), 37°C, 260nm

10µL (0.025mg/mL, DMSO / acetic acid = 100 / 1)