


Unison UK-Phenyl Cadenza CX-C18  
 Unison UK-C8 Cadenza CW-C18  
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Gemfibrozil (Hypolipidemic agent)

ゲムフィブロジル (脂質降下薬)

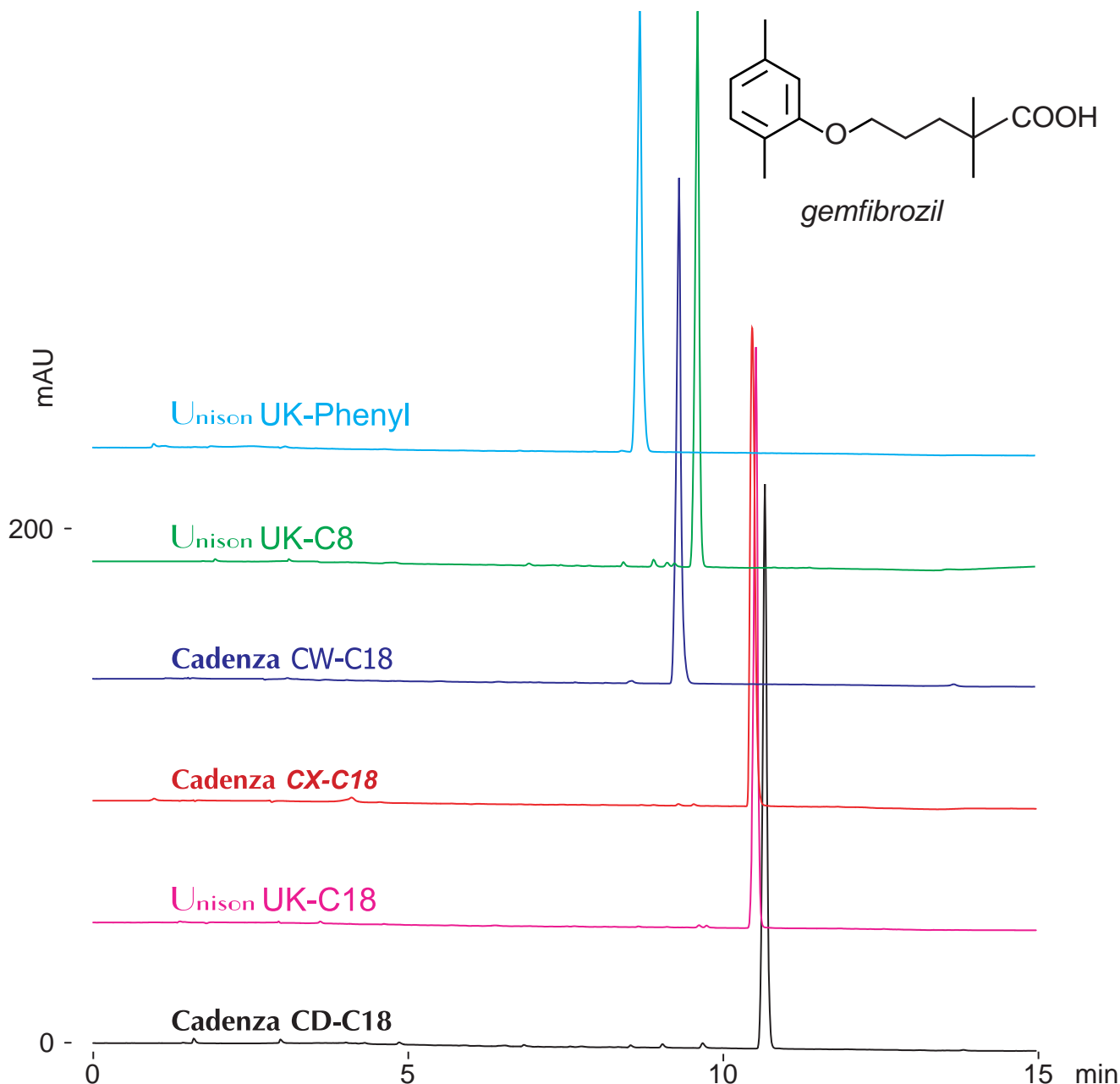

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-12min), 100%B (12-15min)

1.0mL/min (9MPa), 37°C, 280nm

10µL (0.5mg/mL, 100mM ammonium formate / acetonitrile = 70 / 30)