


Cadenza CL-C18  
 Unison UK-C18  
 Cadenza CD-C18

250 x 4.6 mm

HPLC COLUMN APPLICATION

## Cochineal dye (Carminic acid)

コチニール色素 (カルミン酸)

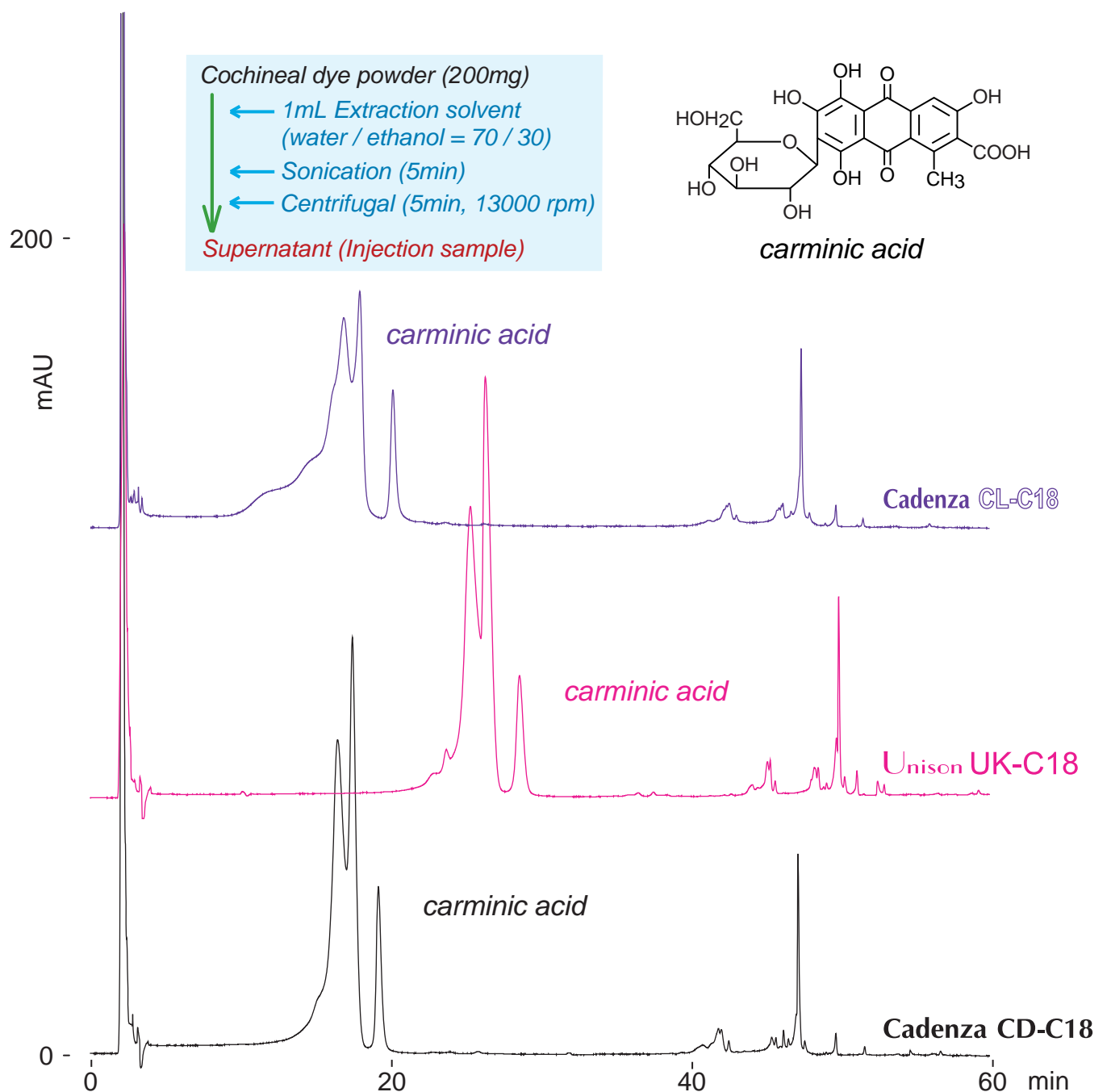

250 x 4.6 mm

A: water / acetic acid = 100 / 0.5

B: water / acetonitrile / acetic acid = 10 / 90 / 0.5

13%B (0-30min), 13-40%B (30-60min)

1.0mL/min (14MPa), 260nm

10µL (extract, water / ethanol = 70 / 30)