



Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Lovastatin (Hypolipidemic agent)

ロバスタチン (高脂血症治療薬)

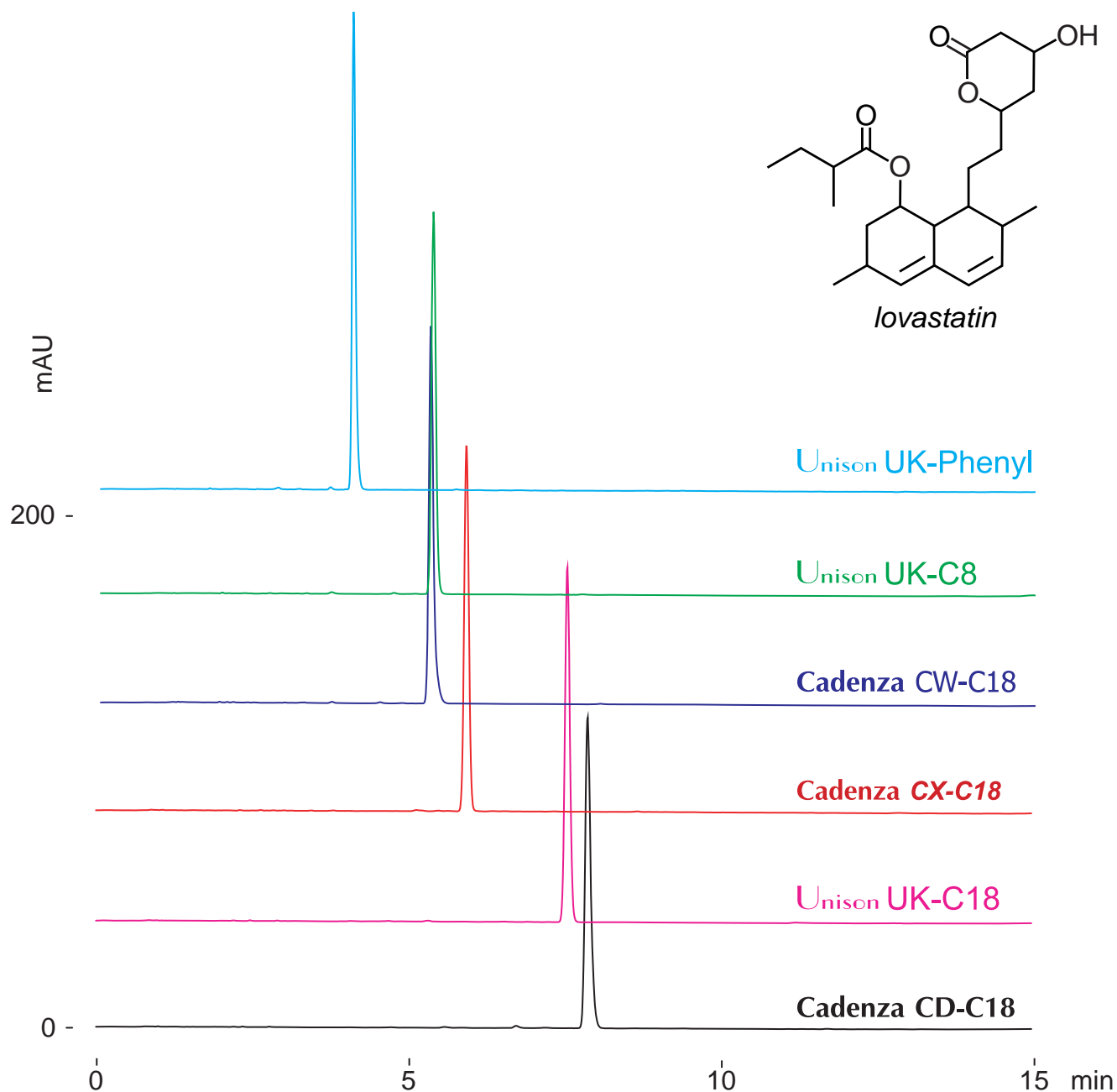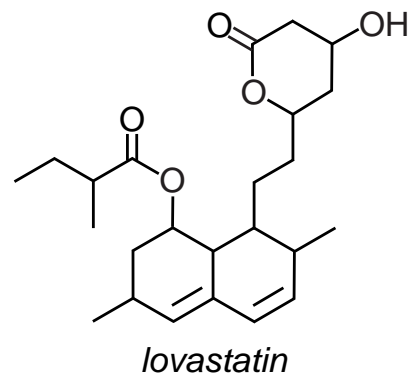

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

50-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 250nm

10µL (0.1mg/mL, 100mM ammonium formate / acetonitrile = 60 / 40)