


Cadenza CD-C18

100 x 3 mm


HPLC COLUMN APPLICATION

LC-MS: Canrenoic acid (Aldosterone antagonist)

カンレノ酸 (アルドステロン拮抗薬) の LC-MS分析

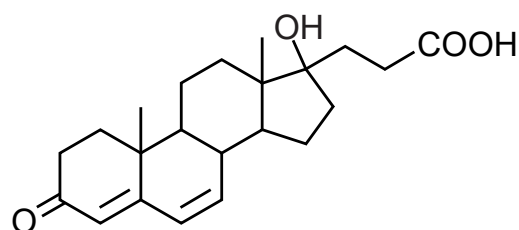

canrenoic acid

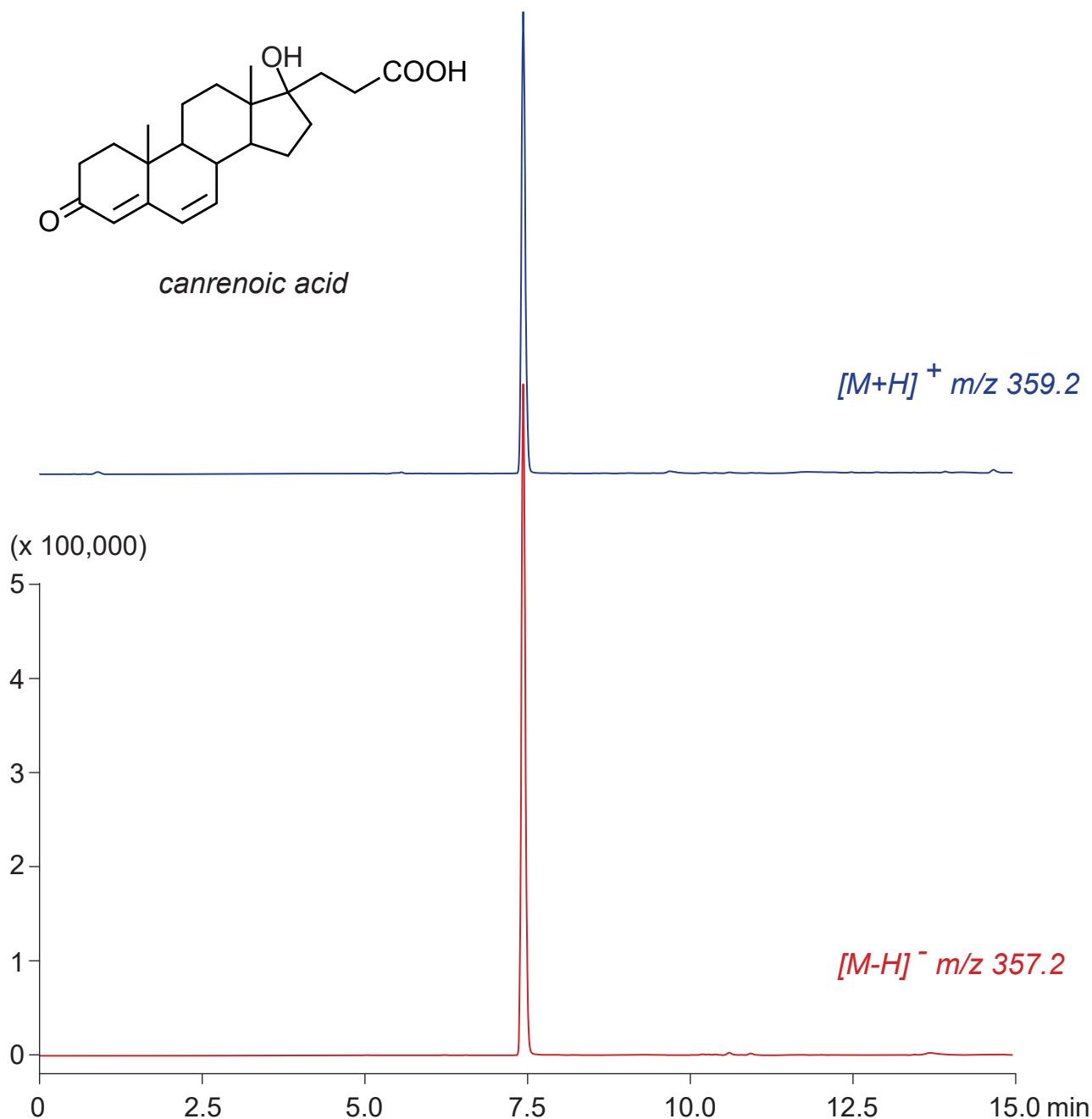

Cadenza CD-C18, 100 x 3 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-100%B (0-12min), 100%B (12-15min)

0.4mL/min (7MPa), 37°C, Single Quad. MS (ESI)

10µL (1µg/mL, 100mM ammonium formate)