


Intrada Organic Acid

100 x 3 mm


HPLC COLUMN APPLICATION

LC-MS: Canrenoic acid (Aldosterone antagonist)

カンレノ酸 (アルドステロン拮抗薬) の LC-MS分析

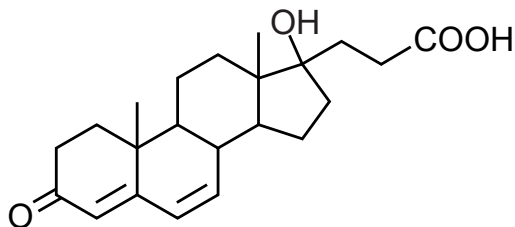

canrenoic acid

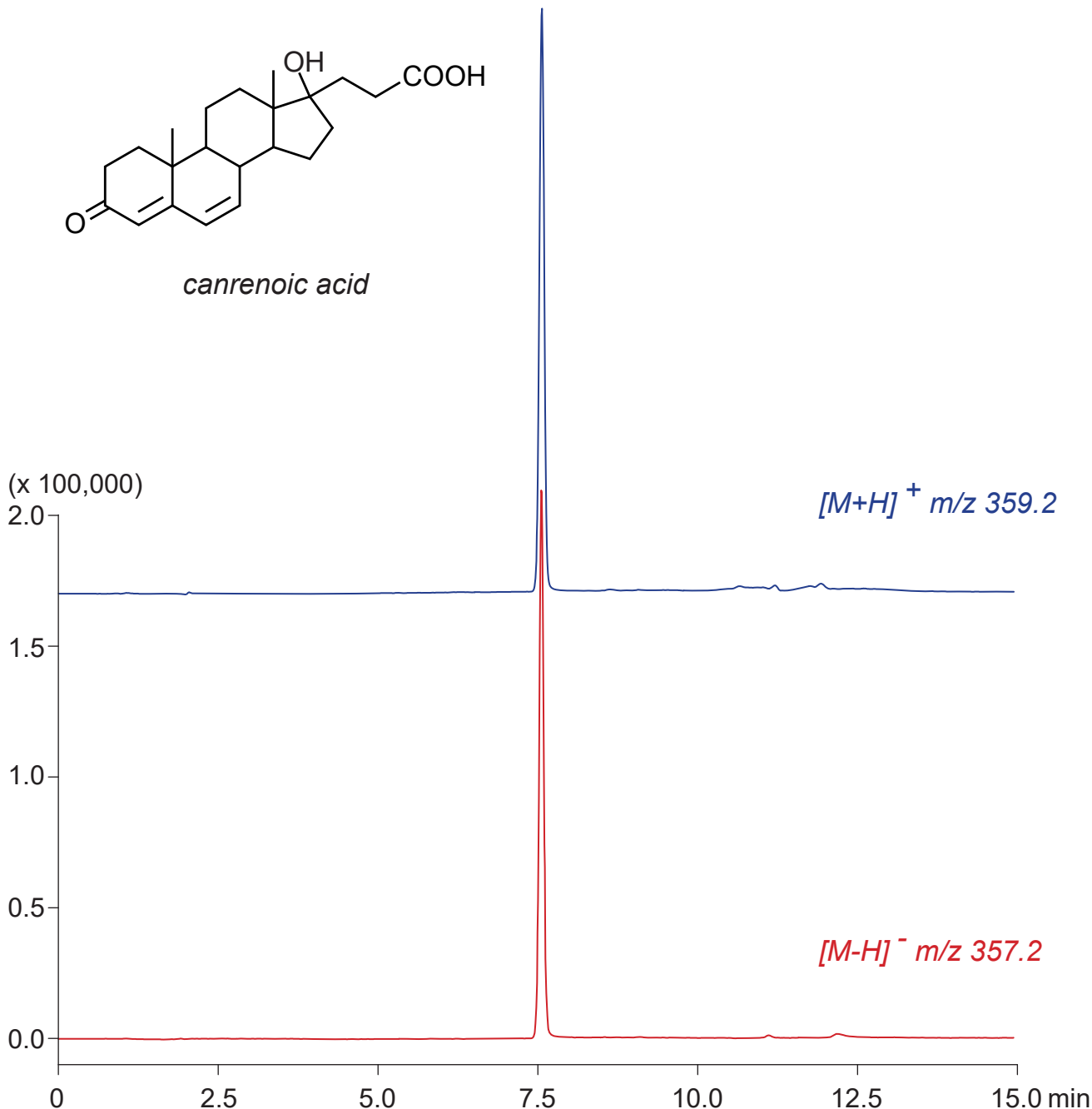

Intrada Organic Acid, 150 x 2 mm

A: water / acetonitrile / formic acid = 80 / 20 / 0.1

B: 100mM ammonium formate / acetonitrile = 40 / 60

0%B (0-1min), 0-100%B (1-12min), 100%B (12-15min)

0.2mL/min (12MPa), 37°C, Single Quad. MS (ESI)

10µL (1µg/mL, 100mM ammonium formate)