

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18

Cadenza CW-C18



Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Triamterene (Potassium-sparing diuretic)

トリアムテレン (カリウム保持性利尿薬)

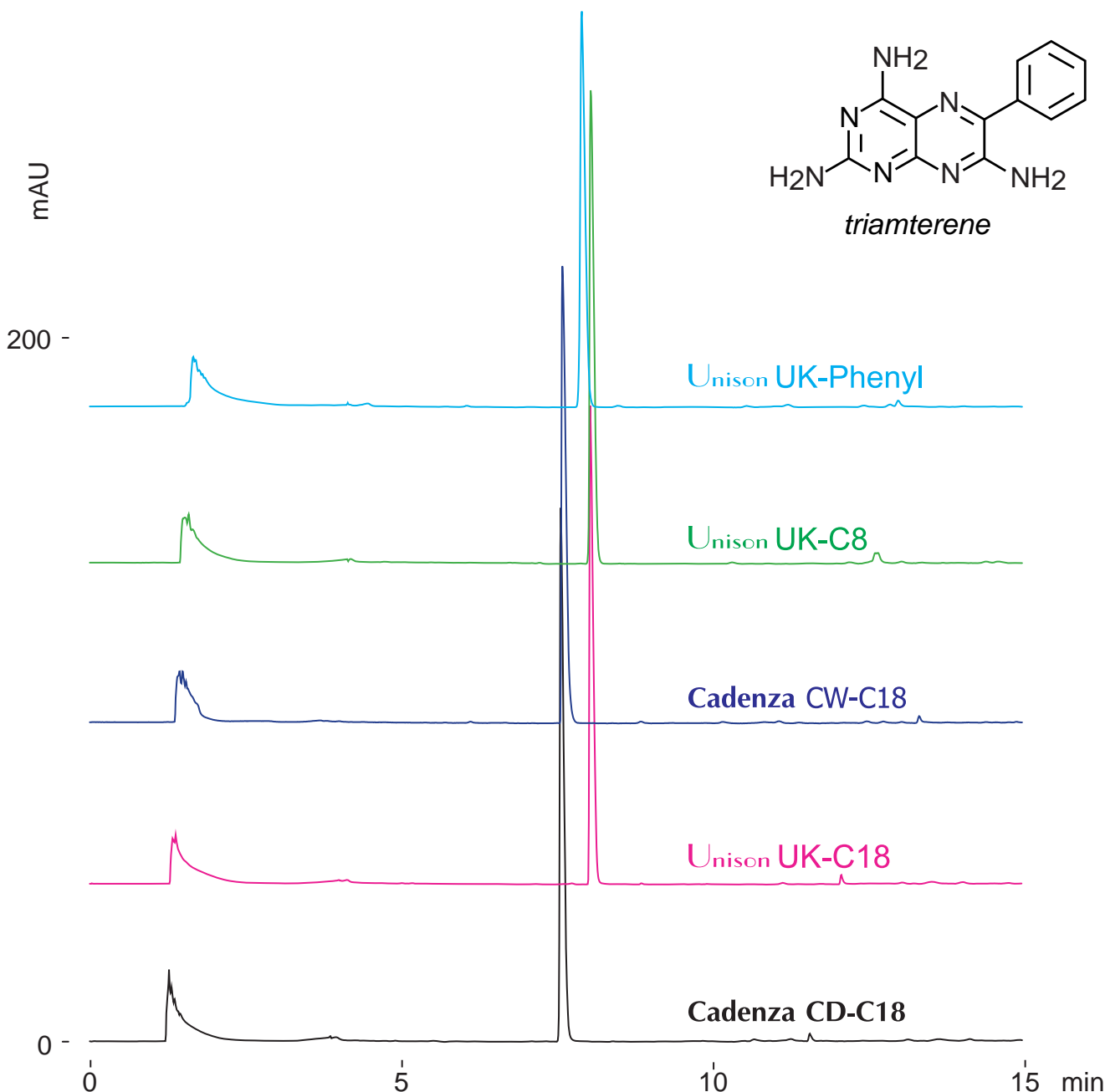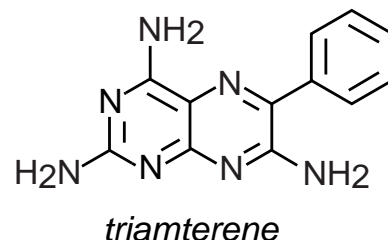

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0%B (0-1min), 0-50%B (1-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.05mg/mL, DMSO / formic acid = 100 / 1)