


Unison UK-Phenyl Cadenza CX-C18  
 Unison UK-C8 Cadenza CW-C18  
 Unison UK-C18 Cadenza CD-C18


100 x 4.6 mm

HPLC COLUMN APPLICATION

## Tris (2,2'-bipyridyl) ruthenium (II) (Dye for dye sensitized solar cell) トリス ( 2,2' -ビピリジル ) ルテニウム ( II ) ( 色素増感太陽電池用色素 )

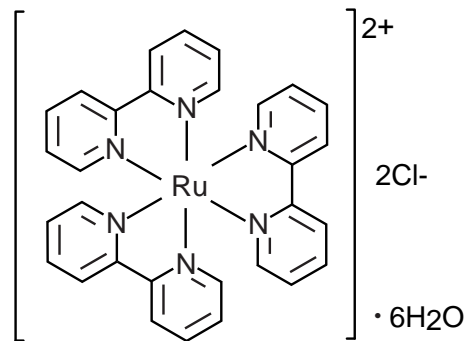

tris (2,2'-bipyridyl) ruthenium (II)  
 chloride hexahydrate

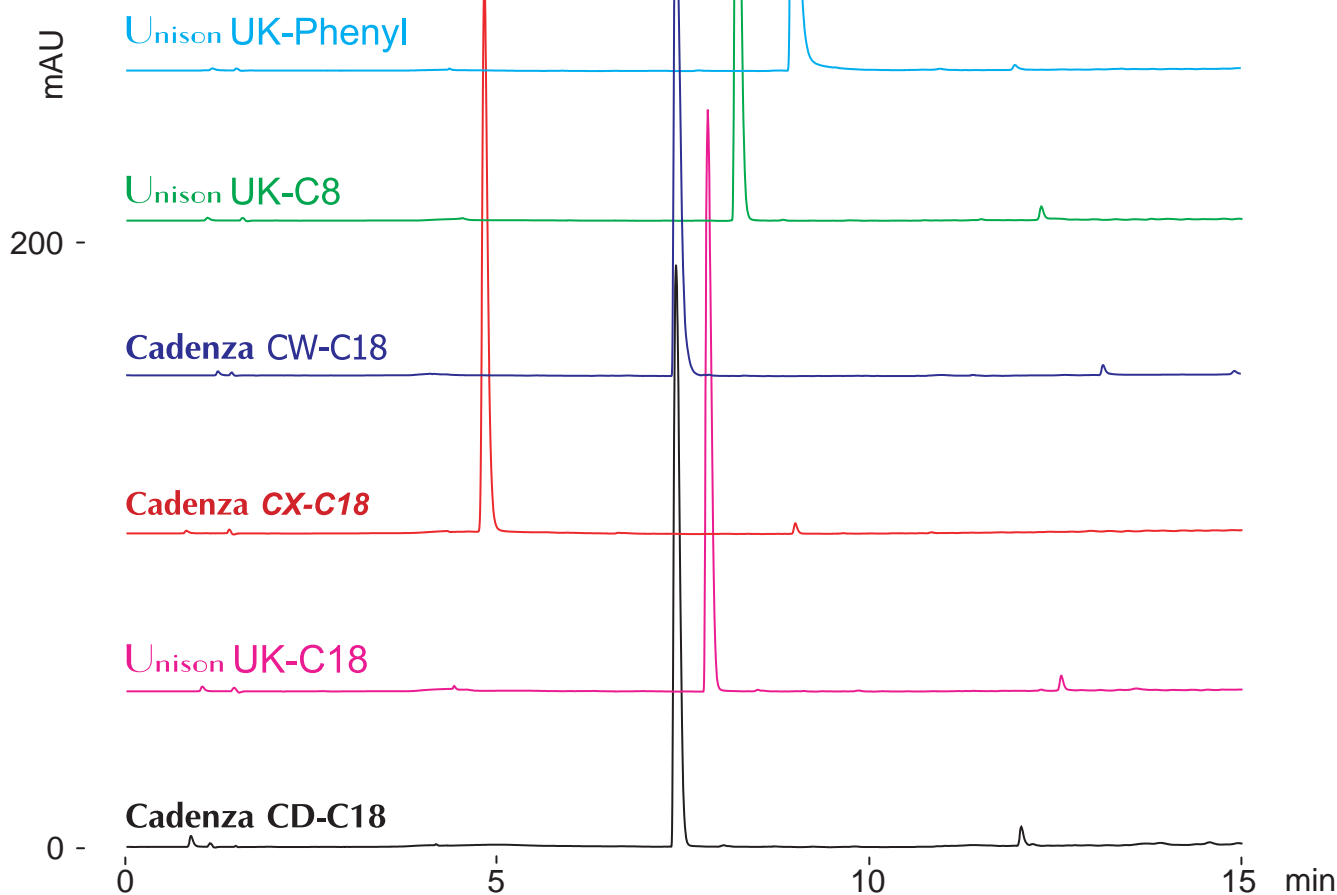

100 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: water / acetonitrile / formic acid = 10 / 90 / 0.5

0%B (0-1min), 0-50%B (1-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.15mg/mL, water / formic acid = 100 / 1)