


Unison UK-Phenyl
 Unison UK-C8
 Unison UK-C18
 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Enzalutamide (Antineoplastic) and synthetic intermediate

エンザルタミド（抗悪性腫瘍薬）と合成中間体

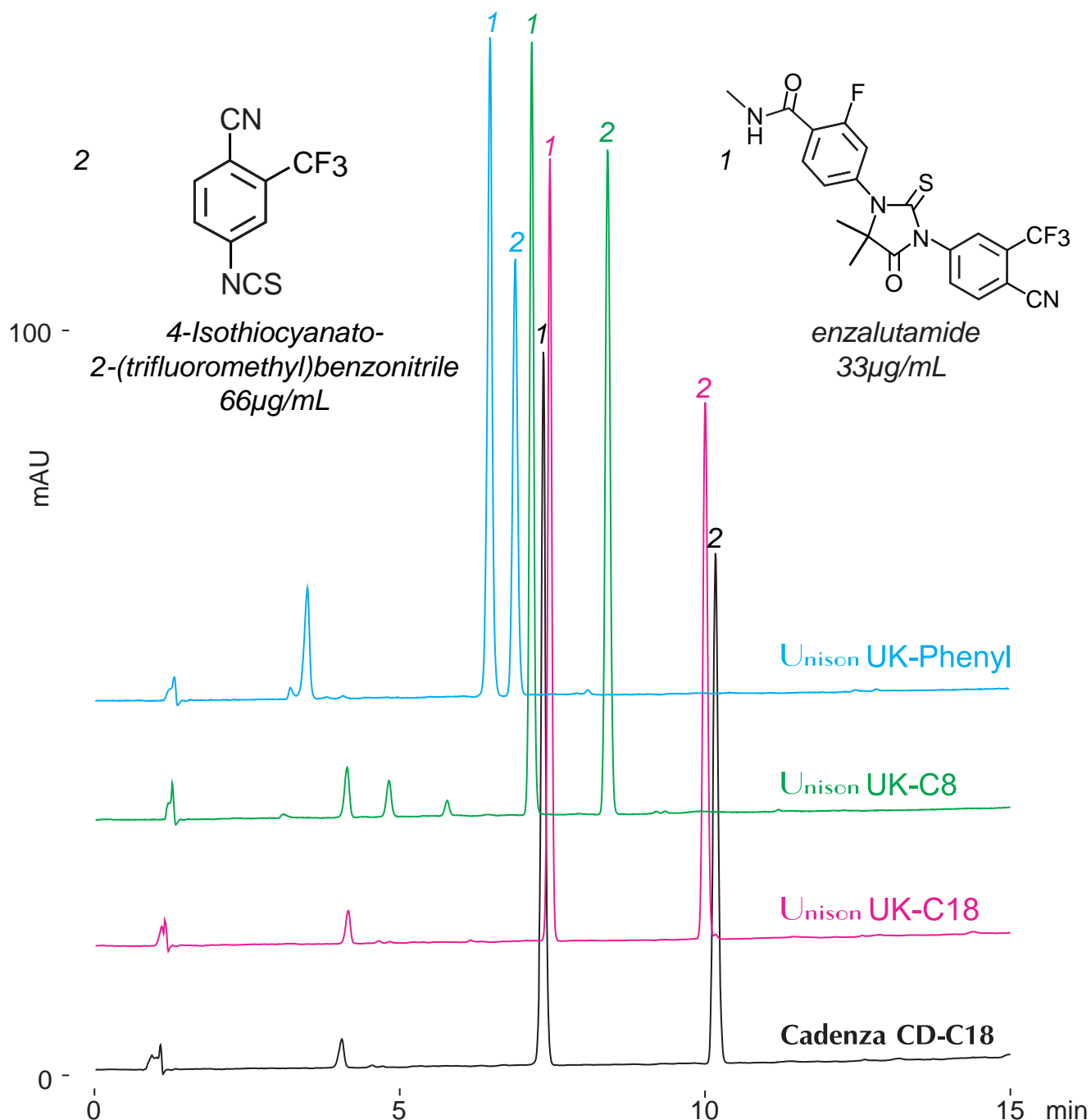

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

30-100%B (0-15min)

1.0mL/min (11MPa), 37°C, 260nm

10 μ L (33-66 μ g/mL, acetonitrile / formic acid = 100 / 1)