




Unison UK-C8

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Sclareol, (+)-sclareolide (Sesquiterpenes)

スクラレオール, (+)-スクラレオリド (セスキテルペン類) の LC-MS 分析

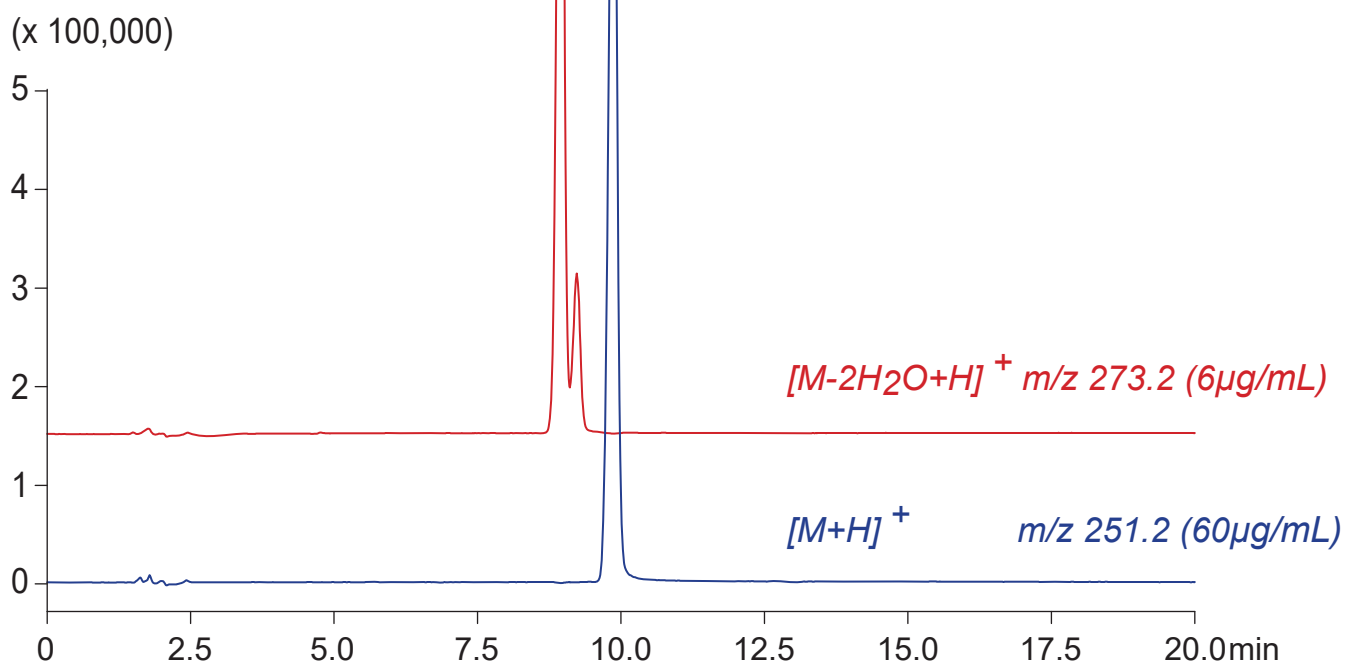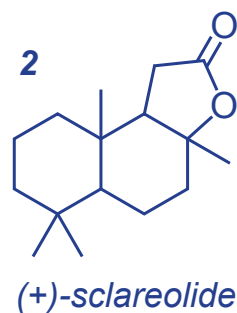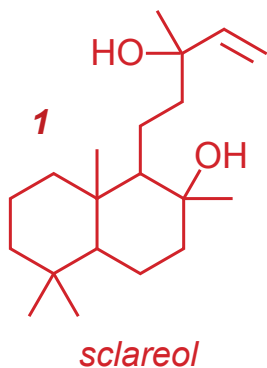

Unison UK-C8, 150 x 3 mm

A: water / acetonitrile / acetic acid = 90 / 10 / 0.1

B: water / acetonitrile / acetic acid = 5 / 95 / 0.1

40-70%B (0-20min)

0.4mL/min (8MPa), 37°C, Single Quad. MS (ESI)

5µL (6-60µg/mL, EtOH / water / AcOH = 50 / 50 / 1)