


Unison UK-C1  
 Unison UK-Phenyl Cadenza CX-C18  
 Unison UK-C8 Cadenza CW-C18  
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

## 1,3-Di-9-carbazolylbenzene (Organic light-emitting diode material)

1,3-ジ-9-カルbazolリルベンゼン (有機エレクトロルミネッセンス材料)

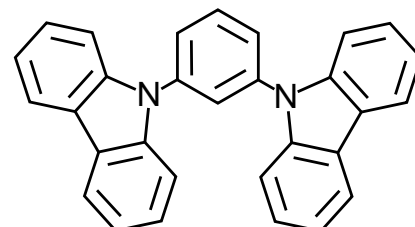

1,3-di-9-carbazolylbenzene

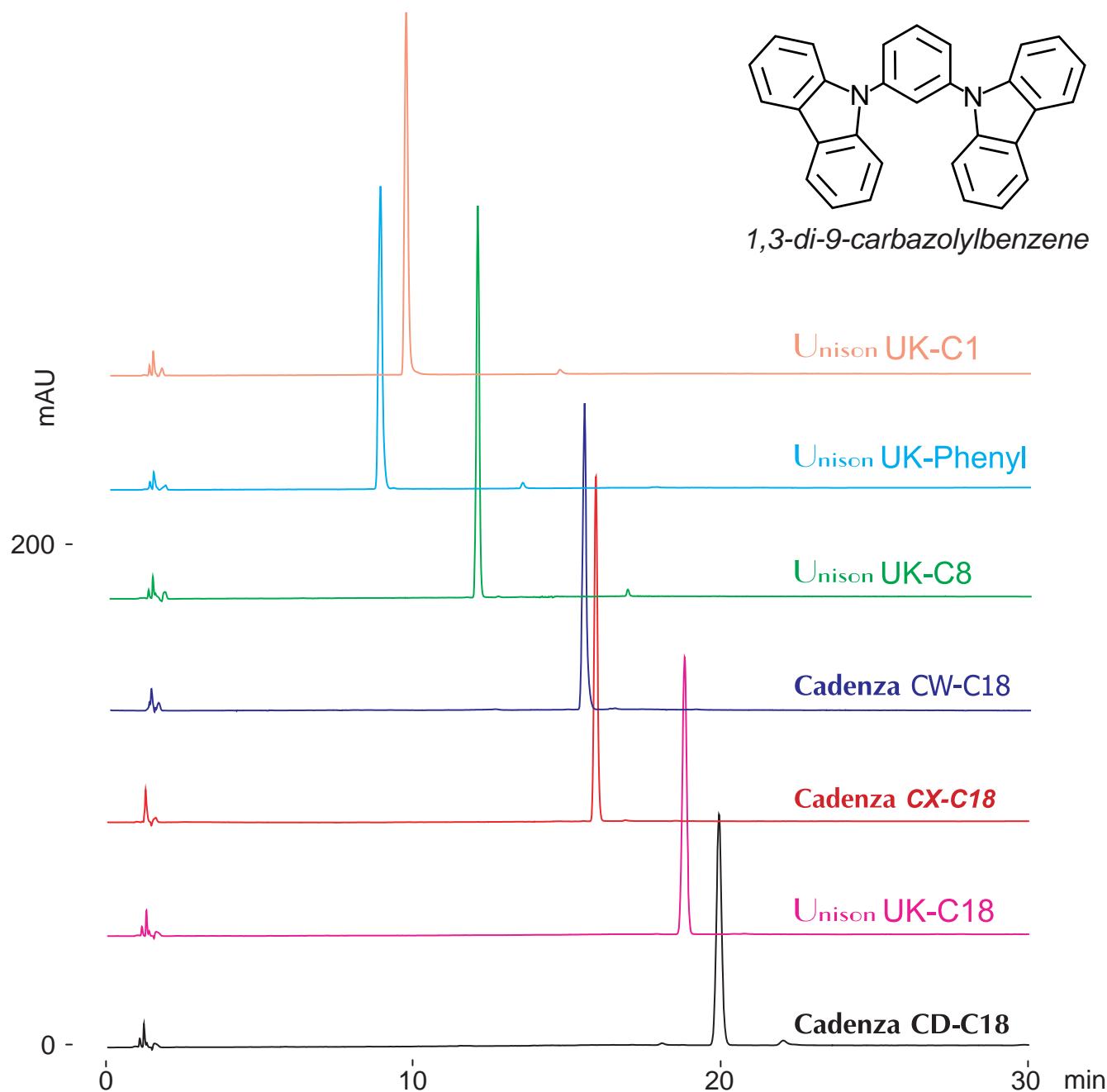

100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-100%B (0-15min), 100%B (15-30min)

1.0mL/min (8MPa), 37°C, 260nm

10µL (50µg/mL, THF / formic acid = 100 / 1)