

Unison UK-C1
 Unison UK-Phenyl
 Unison UK-C8
 Unison UK-C18


Cadenza CX-C18
 Cadenza CW-C18
 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

4,4',4''-Tri-9-carbazolyltriphenylamine (Organic light-emitting diode material)

4,4',4''-トリ-9-カルバゾリルトリフェニルアミン
 (有機エレクトロルミネッセンス材料)

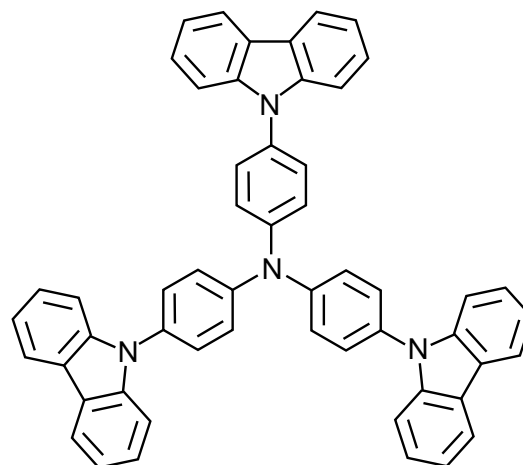

4,4',4''-tri-9-carbazolyltriphenylamine

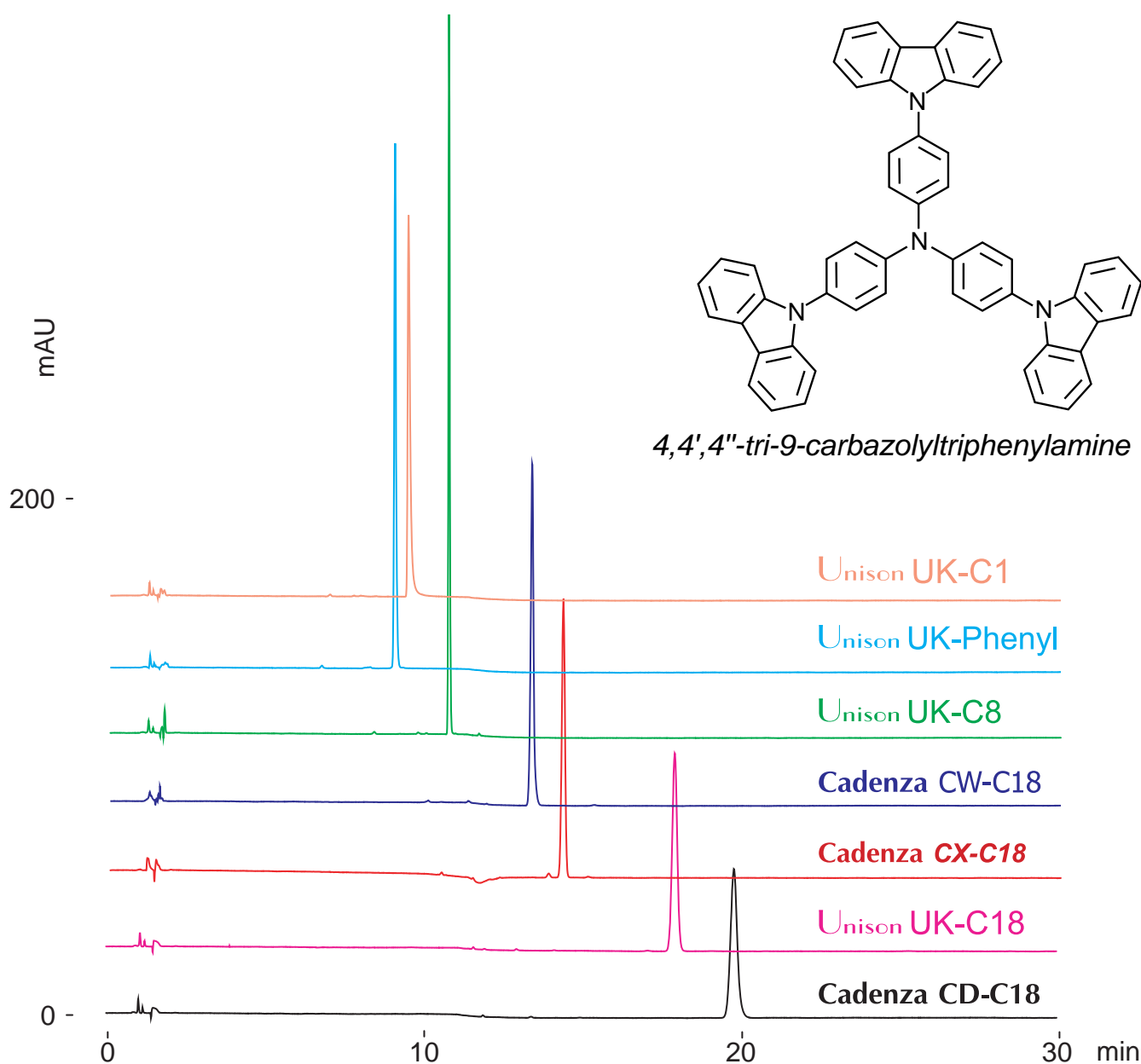

100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-100%B (0-10min), 100%B (10-30min)

1.0mL/min (8MPa), 37°C, 330nm

10µL (0.05mg/mL, THF / formic acid = 100 / 1)