


Unison UK-Phenyl MF 100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: 8-Quinolinol and 8-quinolinolato metal complexes

8キノリノール及び8キノリノラト金属錯体の LC-MS 分析

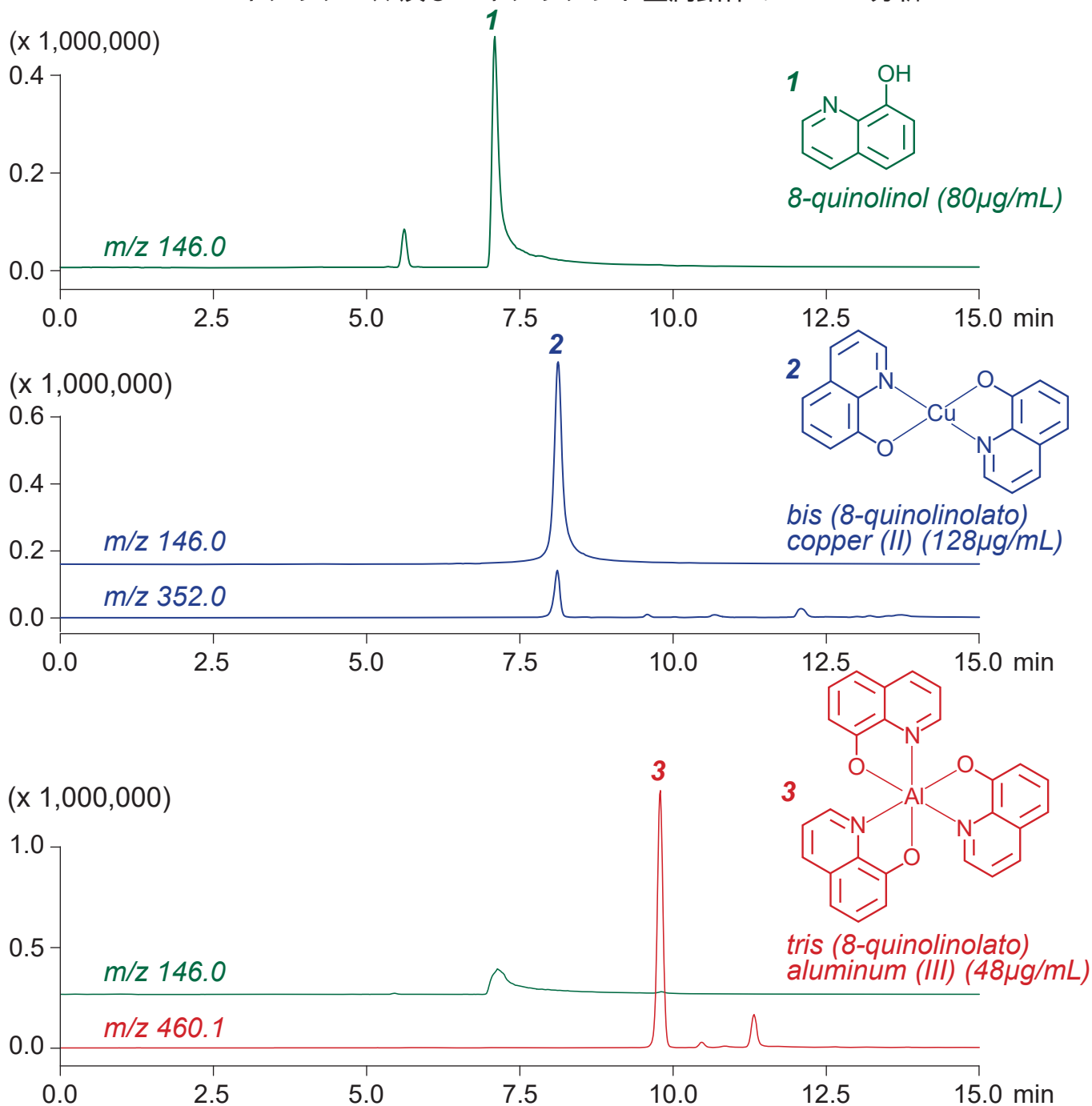

Unison UK-Phenyl MF, 100 x 2 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-15min)

0.2mL/min (9MPa), 37°C, Single Quad. MS (ESI, Positive)

 5µL (48-128µg/mL, 100mM HCOONH₄ / DMSO = 10 / 90)