


Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Coumarin 153 (Organic light-emitting diode material)

クマリン153 (有機エレクトロルミネッセンス材料)

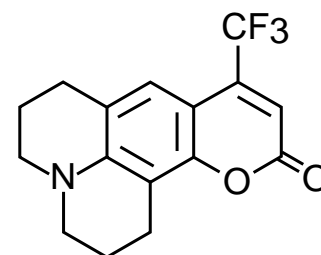

coumarin 153

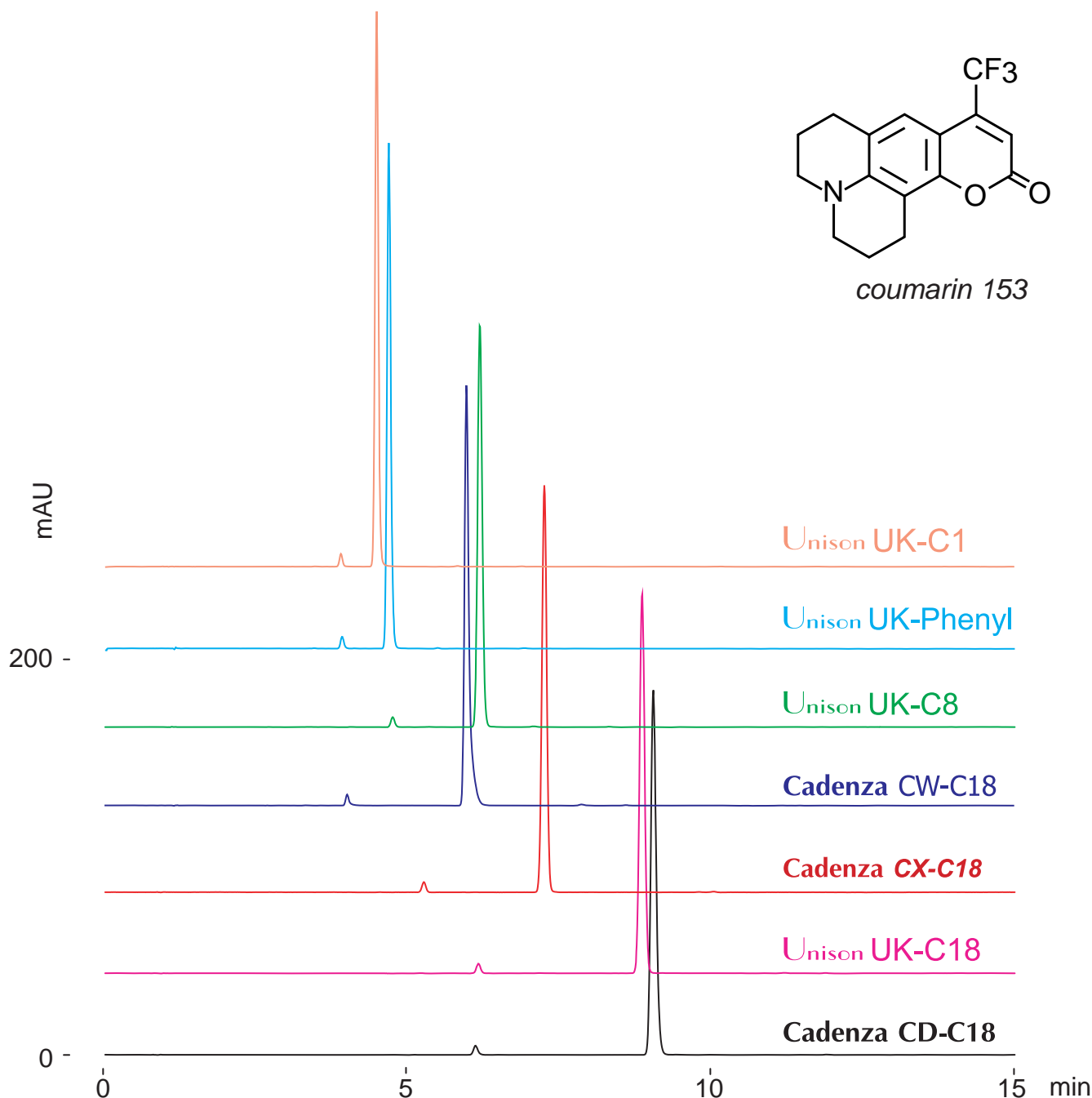

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 50 / 50

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 410nm

10µL (0.1mg/mL, 100mM ammonium formate / acetonitrile = 35 / 65)