


Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Bis [2- (2-benzothiazolyl) phenolato] zinc (II)
 (Organic light-emitting diode material)

ビス [2- (2ベンゾチアゾリル) フェノラト] 亜鉛 (II)
 (有機エレクトロルミネッセンス材料)

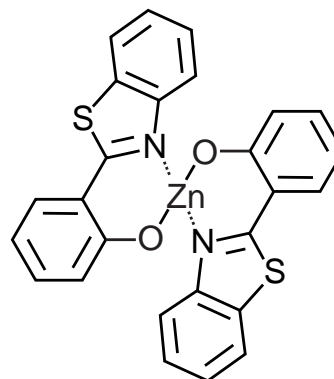

bis [2- (2-benzothiazolyl) phenolato] zinc (II)

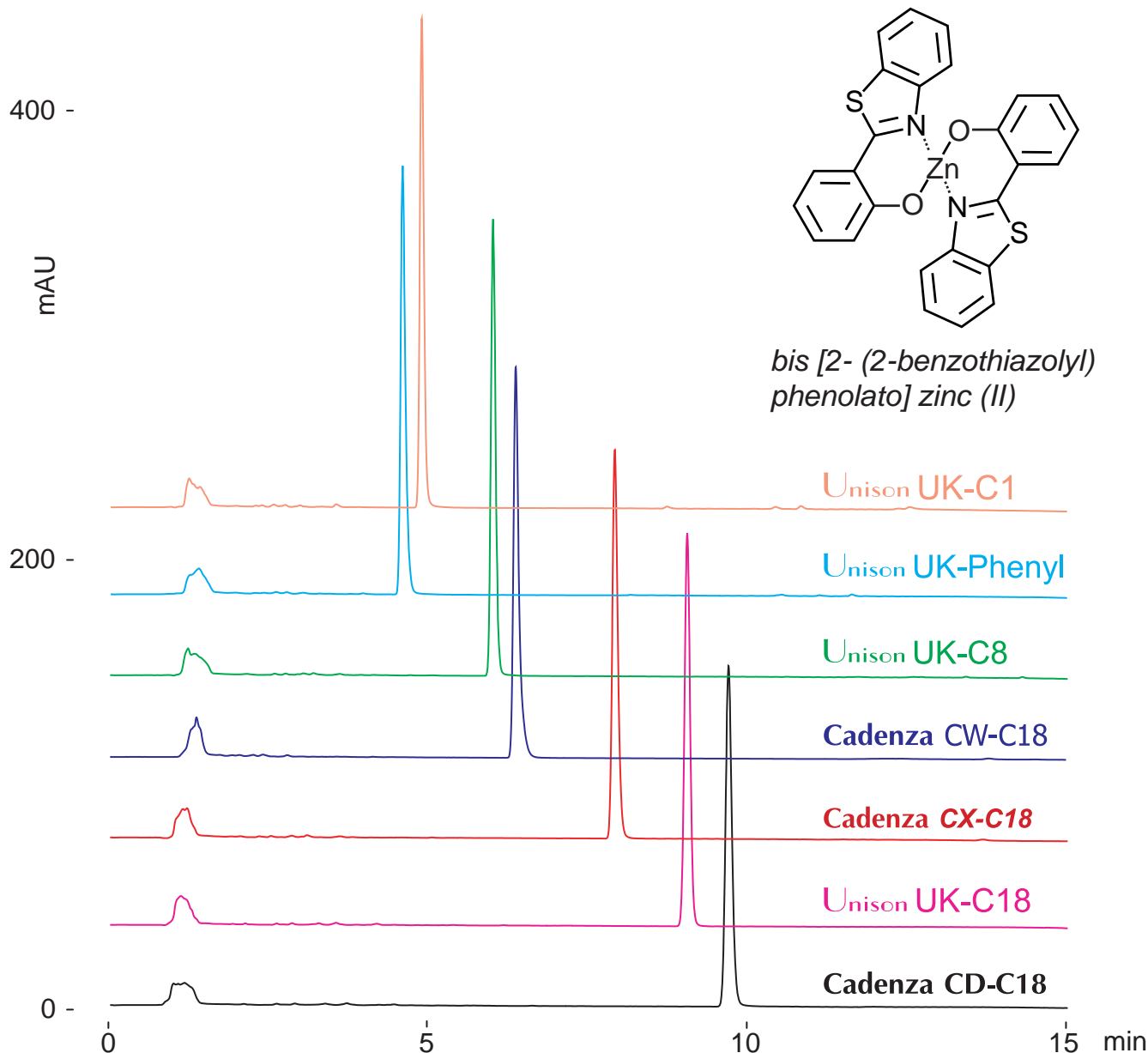

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 50 / 50

B: acetonitrile

0-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 260nm

10µL (0.2mg/mL, 100mM ammonium formate / DMSO = 10 / 90)