



Unison UK-C1  
 Unison UK-Phenyl Cadenza CX-C18  
 Unison UK-C8 Cadenza CW-C18  
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Quinacridone (Organic light-emitting diode material)

キナクリドン (有機エレクトロルミネッセンス材料)

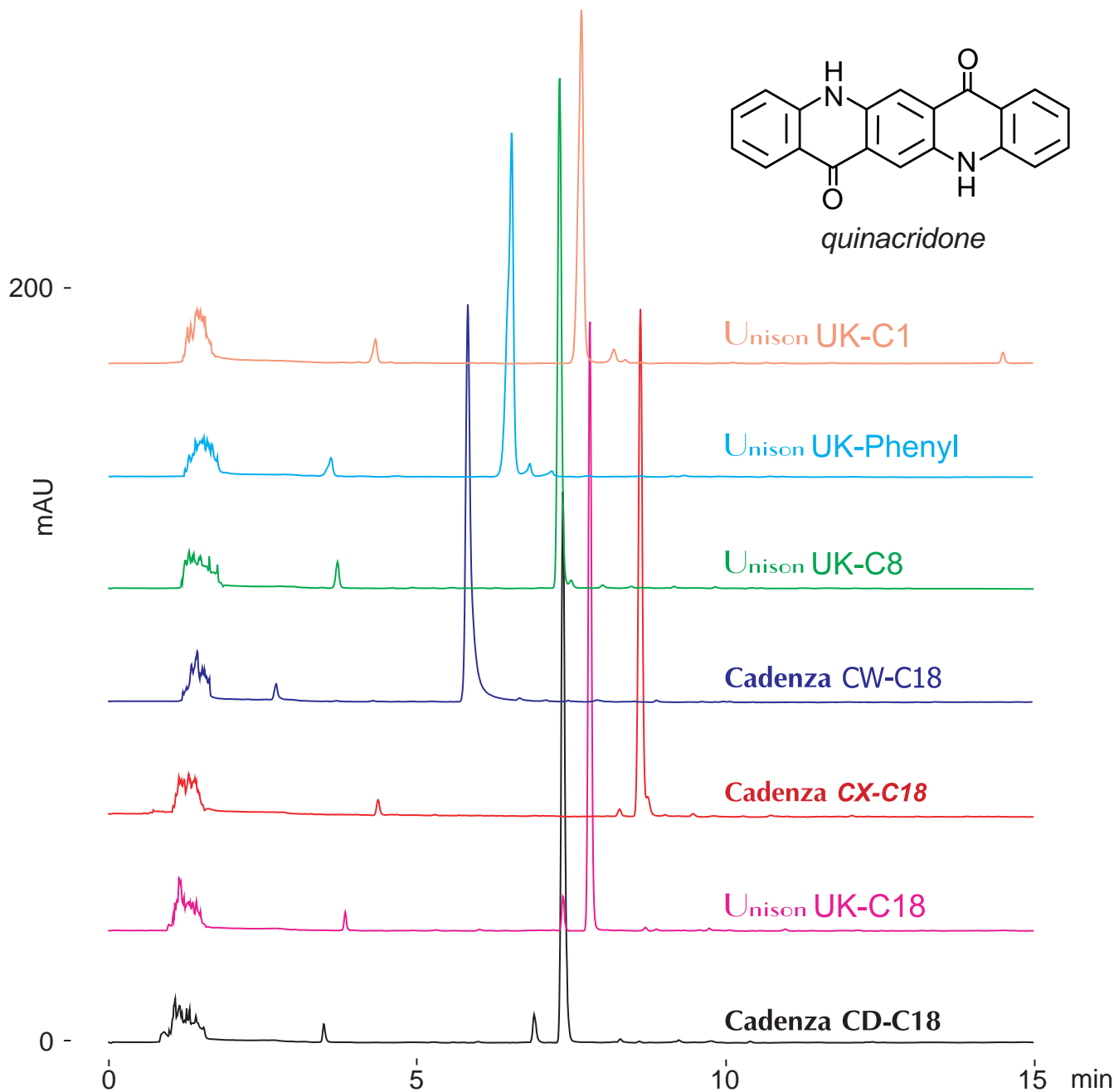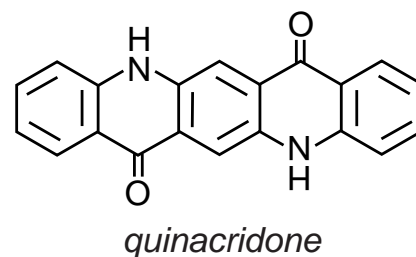

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

10-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 300nm

20µL (12.5µg/mL, DMSO/ formic acid = 100/ 1)