

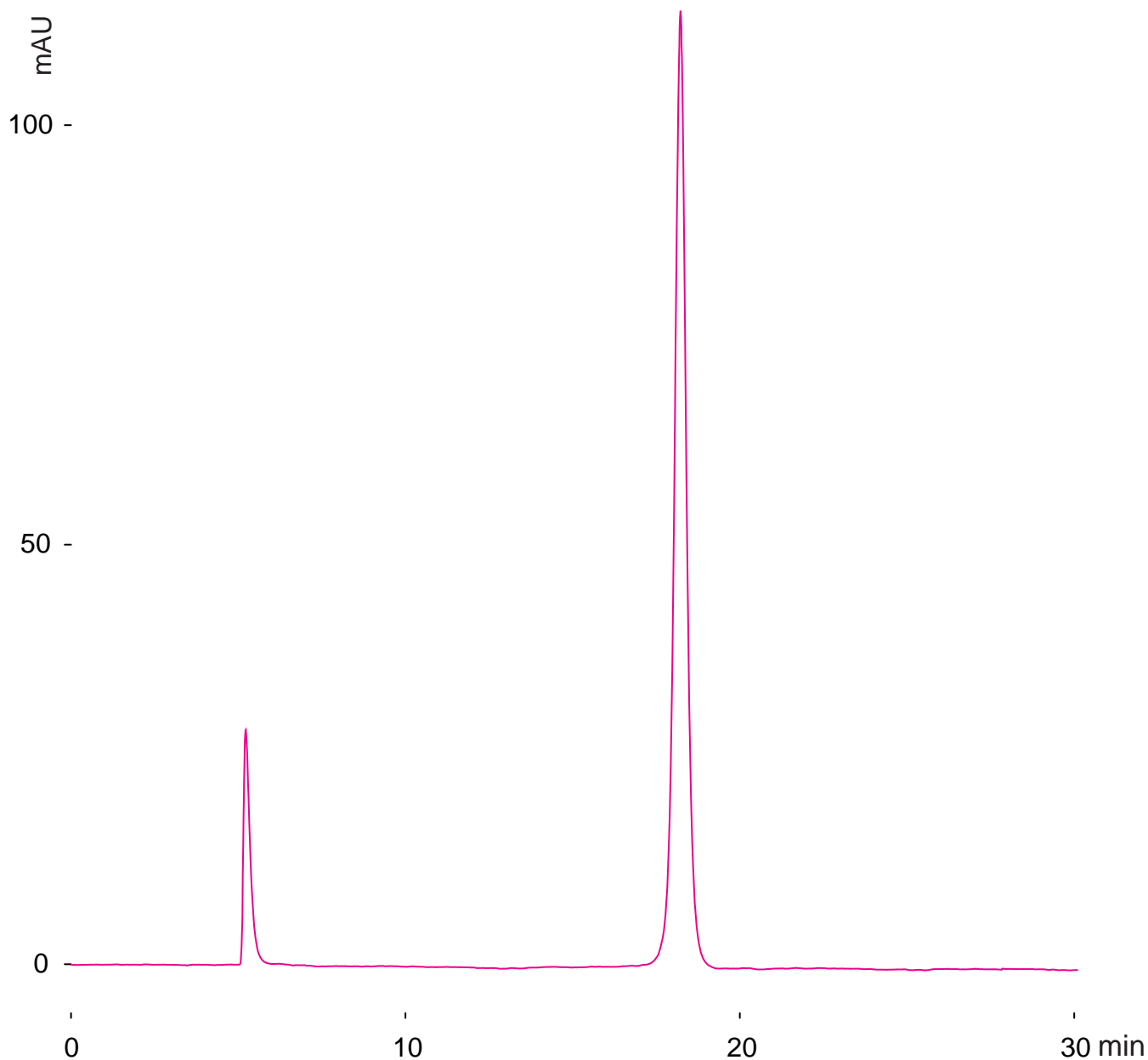
Presto FF-C18

250 x 4.6 mm

HPLC COLUMN APPLICATION

Arabinogalactan from Larch Wood

アラビノガラクトン (カラマツ材由来)

*Presto FF-C18, 250 x 4.6 mm**A: water / formic acid = 100 / 0.1**B: acetonitrile / formic acid = 100 / 0.1**0-20%B (0-30min)**0.3mL/min (22MPa), 37°C**ELSD (spray chamber 70°C, drift chamber 100°C)**10µL (10µg/mL, water)*